

INSTRUCTIONS FOR USE

BREZTRI AEROSPHERE® (brez-TREE)

(budesonide, glycopyrrolate, and formoterol fumarate)

inhalation aerosol, for oral inhalation use

120-inhalation canister

Read this Instructions for Use before you start using BREZTRI AEROSPHERE and each time you get a new refill. There may be new information. This information does not take the place of talking to your doctor about your medical condition or treatment. Your BREZTRI AEROSPHERE inhaler may be different from inhalers you have used before.

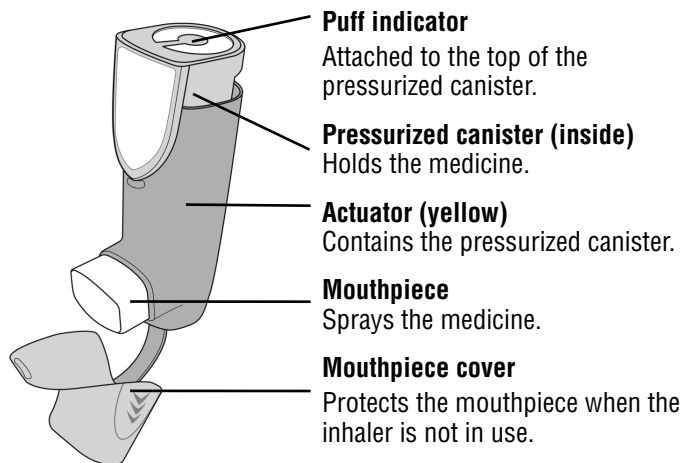
Important information

- **For oral inhalation use only.**
- Prime your BREZTRI AEROSPHERE inhaler before using it for the first time.
- Rinse the yellow actuator 1-time each week.
- Take 2 puffs of medicine in the morning and 2 puffs of medicine in the evening.
- Rinse mouth with water after the 2 puffs to reduce your chance of getting a fungal infection.

Storing your inhaler

- Store your inhaler at room temperature between 68°F to 77°F (20°C to 25°C).
- **Do not store in a humid environment, such as a bathroom.**
- Keep your inhaler and all medicines out of the reach of children.

Parts of your BREZTRI AEROSPHERE inhaler



Reading the puff indicator

The puff indicator will count down by 1 each time you spray a puff of medicine.

Pointer

Points to number of puffs remaining.

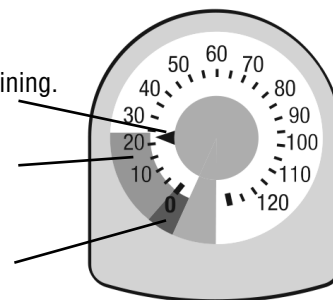
Yellow zone

Order a new inhaler when the pointer is in the yellow zone.

Red zone

Throw away your inhaler when the pointer is at 0 in the red zone.

Do not try to take a puff of medicine when the pointer is at **0** because you will not receive a full dose.




BEFORE FIRST USE – Prime your inhaler 4 times before first use

- Before you use your BREZTRI AEROSPHERE inhaler for the first time, prime it so that you will get the right amount of medicine when you use it.

Prime 1

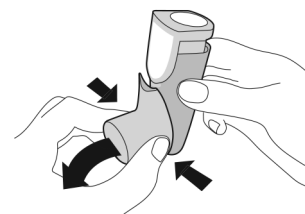
Remove mouthpiece cover by squeezing both sides and pulling out.



Prime 2

Shake the inhaler well and spray **1 test-puff** into the air facing away from you. Repeat Prime 2 for a total of **4 test-puffs**, shaking the BREZTRI AEROSPHERE inhaler before each test-puff.

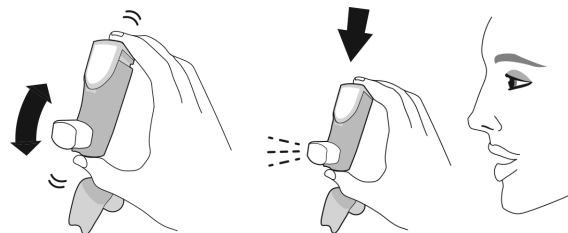**4 total
shake and test-puffs**



Do not skip priming. Extra puffs are provided for priming.

Re-prime your inhaler: To re-prime, spray **2 test-puffs**, shaking before each test-puff.

- after rinsing the yellow actuator
- if dropped
- if not used for more than 7 days

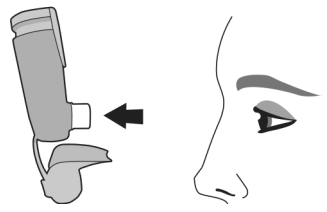
**2 total
shake and test-puffs**

DAILY USE - Using your BREZTRI AEROSPHERE inhaler


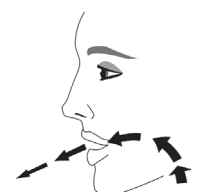

- **Dose: 2 puffs in the morning and 2 puffs in the evening.**
- Rinse mouth with water after the 2 puffs to reduce your chance of getting a fungal infection.

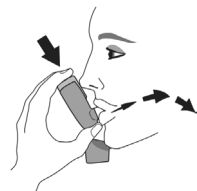
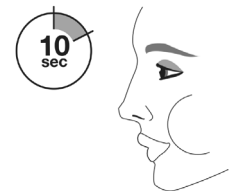
Step 1


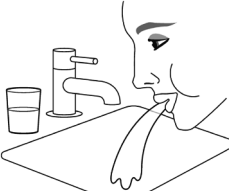
Remove mouthpiece cover. Check the mouthpiece for foreign objects and remove objects before use.



Step 2



<p>2.1 Shake the inhaler well before each puff.</p> 	<p>2.2 Breathe out fully.</p> 	<p>2.3 Place mouthpiece into mouth and close lips around the mouthpiece.</p> 
--	--	---



<p>2.4 Start to breathe in deeply and slowly while spraying 1 puff. Continue breathing in until you cannot anymore.</p> 	<p>2.5 Hold breath for as long as you can, up to 10 seconds.</p> 
---	--



<p style="text-align: center;">Step 3</p> <p style="text-align: center;">Repeat all parts in Step 2 for the second puff.</p>	<p style="text-align: center;">Step 4</p> <p>Put mouthpiece cover back on.</p> 	<p style="text-align: center;">Step 5</p> <p>Rinse your mouth with water. Spit out the water. Do not swallow the water.</p> 
--	---	--


WEEKLY RINSE – Rinse your actuator 1-time each week

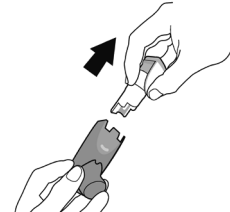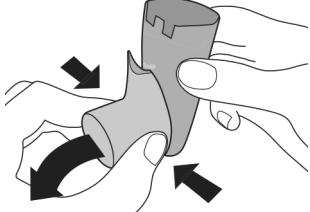- Rinse the yellow actuator 1-time each week so that the medicine does not build up and block the spray through the mouthpiece.
- Keep canister dry.
- Re-prime after rinsing.

<p style="text-align: center;">Rinse 1</p> <p>Remove canister and set aside. Keep canister dry.</p> 	<p style="text-align: center;">Rinse 2</p> <p>Remove the mouthpiece cover.</p> 
--	---
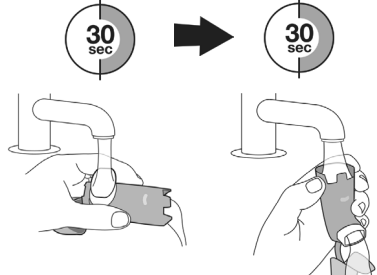
<p style="text-align: center;">Rinse 3</p> <p>Run warm water through the mouthpiece for 30 seconds and then through the top of the actuator for 30 seconds. Rinse for 60 seconds in total.</p> 	<p style="text-align: center;">Rinse 4</p> <p>Shake off as much water as you can.</p>  <p>Do not dry the actuator with a towel or tissue.</p>
---	--
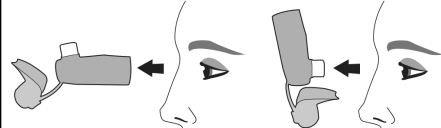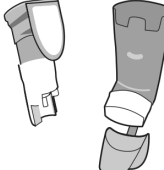
<p style="text-align: center;">Rinse 5</p> <p>Look into the actuator and mouthpiece for medicine build-up. If there is any build-up, repeat Rinse 3, 4, and 5.</p> 	<p style="text-align: center;">Rinse 6</p> <p>Let the actuator air-dry, such as overnight. Do not put the canister back into the actuator if it is still wet.</p> 
--	---

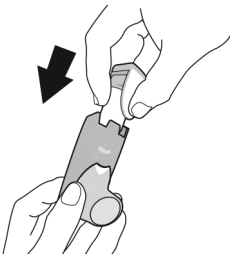Rinse 7	Rinse 8
<p>When completely dry, replace the mouthpiece cover first and then gently press the canister down into the actuator.</p> 	<p>Re-prime the inhaler by spraying 2 test-puffs, shaking before each test-puff.</p> <div style="border: 1px solid black; border-radius: 15px; padding: 5px; text-align: center; width: fit-content; margin: 10px auto;"> <p>2 total shake and test-puffs</p> </div>

Throwing away your BREZTRI AEROSPHERE inhaler

Throw away your inhaler in household trash when:

- puff indicator shows **0**
- or
- **3 months** after your inhaler has been removed from the foil pouch.

Do not reuse or use the actuator with medicine canisters from other inhalers.

Do not puncture or throw the canister into a fire or incinerator.

Ordering a new BREZTRI AEROSPHERE inhaler

- Order a new BREZTRI AEROSPHERE inhaler when the pointer on the puff indicator is in the yellow zone.

BREZTRI AEROSPHERE is a registered trademark of the AstraZeneca group of companies.

© AstraZeneca 2026

Distributed by:
AstraZeneca Pharmaceuticals LP, Wilmington, DE 19850

Revised: April 2026

This Instructions for Use has been approved by the U.S. Food and Drug Administration.

INSTRUCTIONS FOR USE

BREZTRI AEROSPHERE® (brez-TREE) (budesonide, glycopyrrolate, and formoterol fumarate) inhalation aerosol, for oral inhalation use 28-inhalation canister

Read this Instructions for Use before you start using BREZTRI AEROSPHERE and each time you get a new refill. There may be new information. This information does not take the place of talking to your doctor about your medical condition or treatment. Your BREZTRI AEROSPHERE inhaler may be different from inhalers you have used before.

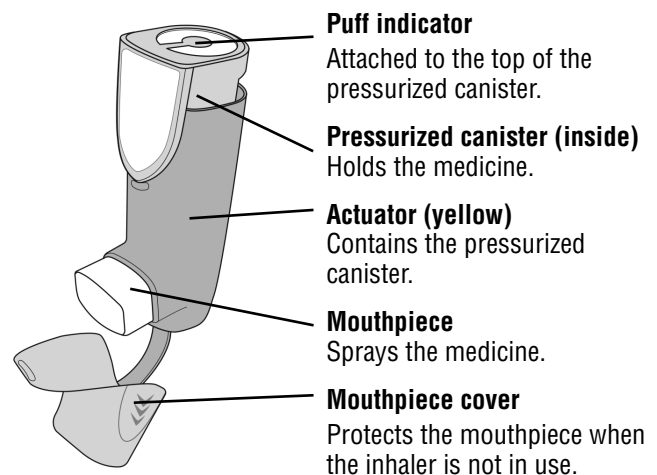
Important information

- **For oral inhalation use only.**
- Prime your BREZTRI AEROSPHERE inhaler before using it for the first time.
- Take 2 puffs of medicine in the morning and 2 puffs of medicine in the evening.
- Rinse mouth with water after the 2 puffs to reduce your chance of getting a fungal infection.

Storing your inhaler

- Store your inhaler at room temperature between 68°F to 77°F (20°C to 25°C).
- **Do not store in a humid environment, such as a bathroom.**
- Keep your inhaler and all medicines out of the reach of children.

Parts of your BREZTRI AEROSPHERE inhaler



Reading the puff indicator

The puff indicator will count down by 1 each time you spray a puff of medicine.

Pointer

Points to number of puffs remaining.

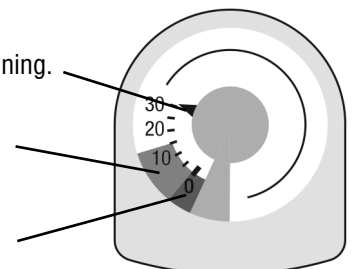
Yellow zone

Order a new inhaler when the pointer is in the yellow zone.

Red zone

Throw away your inhaler when the pointer is at 0 in the red zone.

Do not try to take a puff of medicine when the pointer is at **0** because you will not receive a full dose.




BEFORE FIRST USE – Prime your inhaler 4 times before first use

- Before you use your BREZTRI AEROSPHERE inhaler for the first time, prime it so that you will get the right amount of medicine when you use it.

Prime 1

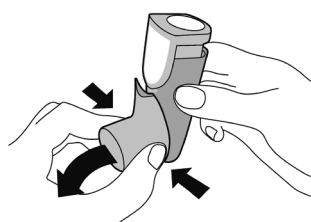
Remove mouthpiece cover by squeezing both sides and pulling out.



Prime 2

Shake the inhaler well and spray **1 test-puff** into the air facing away from you. Repeat Prime 2 for a total of **4 test-puffs**, shaking the BREZTRI AEROSPHERE inhaler before each test-puff.

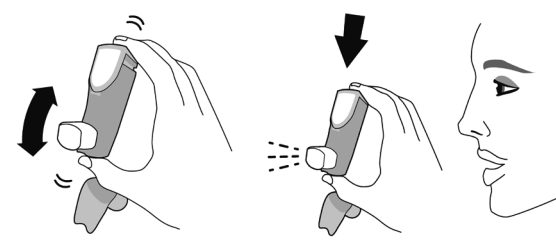4 total shake and test-puffs



Do not skip priming. Extra puffs are provided for priming.

Re-prime your inhaler:

- if dropped
- if not used for more than 7 days

To re-prime, spray **2 test-puffs**, shaking before each test-puff.

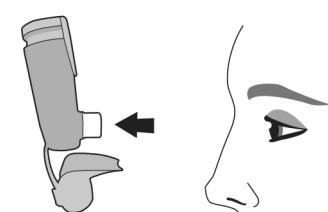
2 total shake and test-puffs

DAILY USE - Using your BREZTRI AEROSPHERE inhaler

- **Dose: 2 puffs in the morning and 2 puffs in the evening.**
- Rinse mouth with water after the 2 puffs to reduce your chance of getting a fungal infection.


Step 1

Remove mouthpiece cover. Check the mouthpiece for foreign objects and remove objects before use.




Step 2

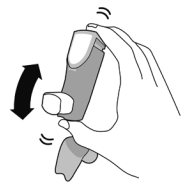
2.1 Shake the inhaler well before each puff.



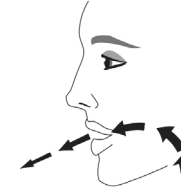
2.2 Breathe out fully.



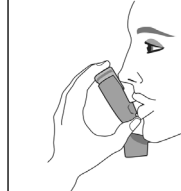
2.3 Place mouthpiece into mouth and close lips around the mouthpiece.





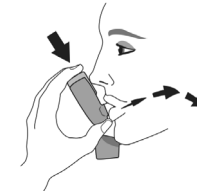2.4 Start to breathe in deeply and slowly while spraying **1 puff**. Continue breathing in until you cannot anymore.



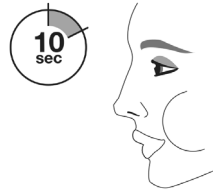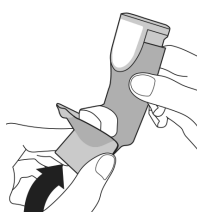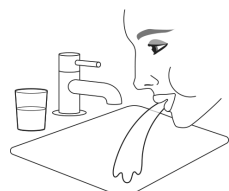2.5 Hold breath for as long as you can, up to **10 seconds**.



Step 3	Step 4	Step 5
<p>Repeat all parts in Step 2 for the second puff.</p>	<p>Put mouthpiece cover back on.</p> 	<p>Rinse your mouth with water. Spit out the water. Do not swallow the water.</p> 

Throwing away your BREZTRI AEROSPHERE inhaler

Throw away your inhaler in household trash when:

- puff indicator shows **0**
- or
- **3 weeks** after your inhaler has been removed from the foil pouch.

Do not reuse or use the actuator with medicine canisters from other inhalers.

Do not puncture or throw the canister into a fire or incinerator.

Ordering a new BREZTRI AEROSPHERE inhaler

- Order a new BREZTRI AEROSPHERE inhaler when the pointer on the puff indicator is in the yellow zone.

BREZTRI AEROSPHERE is a registered trademark of the AstraZeneca group of companies.
©AstraZeneca 2026

Distributed by:
AstraZeneca Pharmaceuticals LP, Wilmington, DE 19850

Revised: April 2026

This Instructions for Use has been approved by the U.S. Food and Drug Administration.